



THE SOLUTION OF CHOICE
*for the treatment of
Dysfunctional Uterine Bleeding*

PNN MEDICAL SA

***Cooperating with healthcare
professionals to improve patient's quality of life***

Pnn Medical SA is committed to improve the quality of life for women by offering minimally invasive therapies.

Using best practice methods in research, development, production, distribution, and product training we are providing world wide access to Gold Standard

medical equipment with a particular focus on product quality, patient safety, and procedure efficacy.

Our goal is to support healthcare professionals optimising clinical results and satisfaction for women treated with our products.



Distributed by



Pnn Medical SA

Av. Riond-Bosson 14 | CH – 1110 Morges
Tel. +41 21 804 99 00 | Fax +41 21 804 99 01
info@cavaterm.com | www.cavaterm.com



Cavaterm™ is an endometrial ablation system for the permanent treatment of Dysfunctional Uterine Bleeding. It consists of a single use catheter and a Cavabox.

The combined action of circulating heat and adjustable pressure inside the uterine cavity during 10 minutes will uniformly and irreversibly coagulate


the endometrium and underlying myometrium. Specific features built in the system ensure maximal safety and efficiency.

Your patient will benefit from a safe and minimally invasive solution, which will significantly improve her quality of life with relieve from the bleeding disorders.


More than 90% of the patients treated with **Cavaterm™** worldwide are satisfied and the success rates obtained are maintained over the years.

THE Cavaterm™ TREATMENT


Balanced combination of temperature, pressure and time ensures permanent treatment of heavy menstrual bleeding.

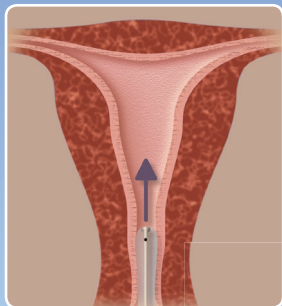

Insertion of a balloon catheter into the uterine cavity

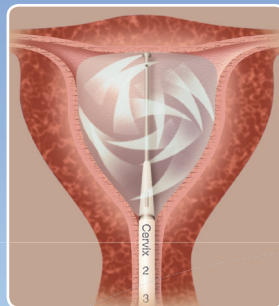

Application of circulating heat and pressure

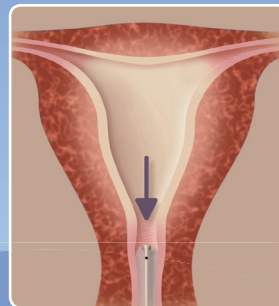

Removal of the catheter

FEATURES

A clever combination of unique features for proven efficacy and safety

- The flexible and adjustable silicone balloon is very conforming and optimally fits the uterine cavity shape and size.
- The small 6 mm catheter diameter allows easy device insertion.
- Cervical insulation protects cervical canal and vagina from heat effects.
- The self-regulating heating elements are safe and powerful heat providers, which ensure maximal protection and efficient heat treatment.
- Vigorous fluid circulation ensures an even heat diffusion inside the uterine cavity for a uniform heat treatment.
- Constant pressure can be easily and precisely maintained at optimal level throughout the treatment with a special Pressure Regulating-Syringe, for efficient heat penetration through the endometrial tissue without risk of overpressure.
- The Cavabox is light-weight, compact and extremely easy to operate, making set-up, use, transport and storage simple, quick and convenient.
- The system is equipped with a memory chip retaining treatment data.

BENEFITS

Improved quality of life through a safe and minimally invasive solution

- Outpatient treatment
- Simple and short
- Permanent solution
- Free of hormones
- High acceptability among patients
- Success and satisfaction rates of more than 90 %
- Swiss quality product

